# ARDESTO™

## ELECTRIC TOOTHBRUSH

Instruction manual

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**ETB-B300BL, ETB-B300W**

Read the instruction manual carefully before using this appliance – electric toothbrush. The instructions below will help to use the appliance effectively, hygienically and safely. Save this manual for future reference.

## **SAFETY MEASURES**

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- This appliance is intended only for the personal hygiene of the human oral cavity.
- A toothbrush is an item of personal hygiene and should be used individually.
- Consult a dentist before using a toothbrush, especially if oral or gum surgery have had a place in the past two months, have undergone dental treatment, have dental problems, or have a pacemaker or other implants.
- When using this appliance for the first time, due to the high-speed vibration of the bristles, some people may experience irritation or numbness of the teeth, and may experience slight bleeding of the gums or mild discomfort. This is a normal reaction. If such symptoms persist after two weeks of regular use, discontinue use and consult a dentist.
- If experience pain during use, stop use immediately and consult a dentist.
- If the appliance has been stored or transported at low temperatures, keep the appliance at room temperature for 2 hours before use.

- The appliance must be powered only from a safe extra-low voltage according to the marking on the appliance.
- Before first use and regularly thereafter, check the appliance and cable for damage. Do not use the appliance and the cable if they are not working properly or if they have been dropped or damaged in any other way. Do not attempt to repair the appliance and accessories by yourself. Contact an authorized service center for diagnosis or repair.
- Do not attempt to change or modify the appliance. Do not insert any objects into the openings of the appliance.
- Use only accessories and attachments supplied with the appliance or those indicated as suitable for use with this appliance or accessories provided by an authorized service center.
- Do not use the appliance without the toothbrush head.
- If a toothpaste contains impurities – peroxide, baking soda or other bicarbonates (often used in whitening toothpastes) thoroughly wash the toothbrush head and toothbrush body after use to avoid damaging of the appliance plastic with residual impurities.
- Do not damage, modify the power cable, and do not stretch, bend or twist the power cable. Do not press heavy objects on the power cable or get it caught between objects.
- Avoid contact of the appliance and power cable with heat sources.
- Do not store the appliance in places with high temperatures, avoid direct sunlight.
- Do not use the appliance while charging.
- Do not charge the appliance in humid or high temperature locations.
- Make sure the charging connector and power cable are never wet or damp.

- Do not immerse the appliance in water or other liquids.
- The appliance is intended for domestic use only as intended in accordance with this manual. The appliance is not intended for use in commercial, industrial or medical environments.
- Do not use the appliance outdoors. Do not use the appliance outdoors or in areas with high oxygen levels or where aerosols (sprays) are used. The appliance is sensitive to chemicals and high oxygen levels and aerosols create dangerous conditions for the electronics and can cause short circuits, damage the appliance, and fire or explosion.
- Do not disassemble the appliance to remove the battery. The appliance has a built-in battery. Do not attempt to remove it from the appliance yourself, contact an authorized service center.
- Leakage of electrolyte from the battery can be caused by overloading the appliance or using it at extremely high temperatures. If battery fluid gets on the skin or eyes, the affected area should be washed immediately with water. Seek medical attention.
- Use the appliance within the permissible ambient temperature range – from 0 to 40°C.

## EQUIPMENT

**Electric toothbrush** (body) – 1 pc

**Toothbrush heads** – 2 pcs

**Travel case** – 1 pc

**USB-A Type-C cable** – 1 pc

**Toothbrush holder** – 1 pc

**Instruction manual (contains a warranty card)** – 1 pc

## APPLIANCE DESCRIPTION


**Power supply:** li-ion battery 3,7 V, 500 mAh, 14500


**USB-A Type-C cable** compatible with any 5 V  $\overline{\text{---}}$  1 A adapter

**Waterproof class:** IPX7

The manufacturer reserves the right to make changes to the technical specifications and design of appliances that do not impair performance without prior notice.


Pictures in the manual are for reference only. Actual appliance may vary.
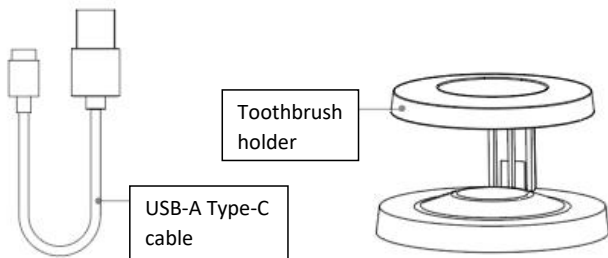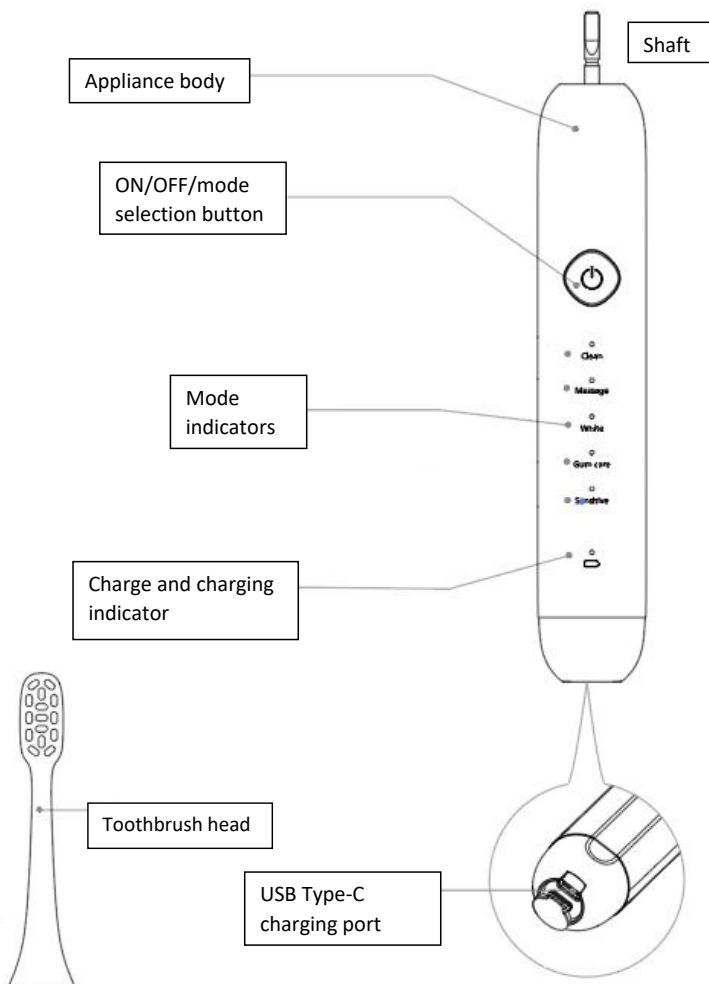
## **APPLIANCE CHARGING**

1. The red flashing charge indicator or the appliance not turning on means that the appliance battery needs recharging.
2. Open the waterproof silicone plug.
3. Connect the USB-A Type-C cable to the port on the appliance and to a power source (to the charger port).
4. When the charge indicator flashes green slowly, charging is in progress.
5. When the charge indicator light is solid green for 30 seconds and then turns off, the battery is fully charged.
6. After charging is complete, unplug the appliance and close the silicone plug.

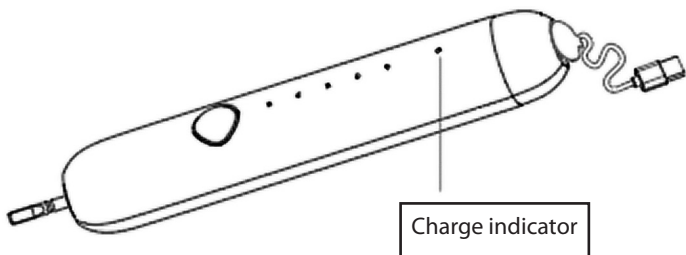
A full charge takes approximately 3 hours. Do not use the appliance while charging.

If the appliance has not been used for a long time and the indicator does not turn on when recharging, wait 6-10 minutes for recovery.


The appliance body may be slightly warm during charging, this is normal. If the body becomes hot, turn off the power immediately and contact an authorized service center.

## APPLIANCE USAGE


1. Place the toothbrush head on the appliance body shaft. Please the toothbrush head lightly to lock in onto the toothbrush body. There is a small gap between the toothbrush head and appliance body.

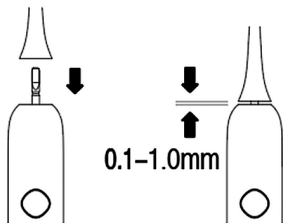
2. Wet the cleaning part of the toothbrush head and apply a small amount of toothpaste.

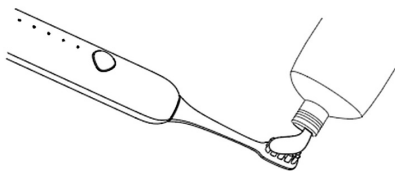3. Hold the toothbrush head at a slight angle to the teeth and press down it touches the gum line.

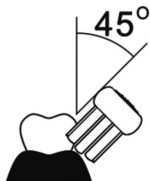4. Press ON button.
5. Briefly press this button again to select the required mode (the corresponding mode indicator will light up):
  - **Clean:** standard setting for effective plaque removal and general oral hygiene.
  - **Massage:** delicate cleaning with a massage effect for the gums, gentle stimulation mode to improve blood circulation.
  - **White (Whitening):** intensive cleaning to remove surface stains from enamel, making teeth whiter.
  - **Gum care:** gentle movements to clean the gum line without excessive pressure, preventing inflammation and bleeding.
  - **Sensitive:** delicate yet thorough cleaning for sensitive teeth and gums.
6. Do not use excessive force when brushing, press the toothbrush head lightly against the teeth moving it gently back and forth.
7. To ensure even brushing divide the mouth into 4 areas. Brush each area for 30 seconds. Start with area 1 (outside of upper teeth), then move to area 2 (outside of lower teeth), then move to area 3 (inside of upper teeth), and finally


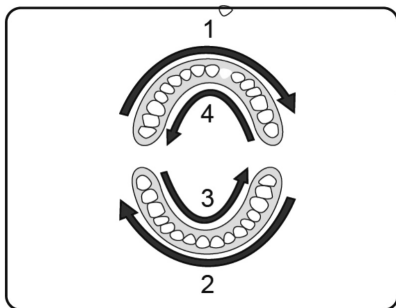
to area 4 (inside of lower teeth). During cleaning, a slight vibration warning is felt every 30 seconds, indicating the need to move on to the next area.

8. The appliance will automatically turn off after 2 minutes of use. To turn off the appliance at any time, briefly press the button until it turns off.

## **CLEANING AND CARE**

- Be sure to disconnect the appliance from the power source before cleaning.
- Clean the appliance after each use.
- Do not immerse the appliance body and cable in water or other liquids.
- Clean the power cable with a dry, soft cloth.

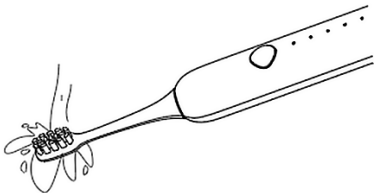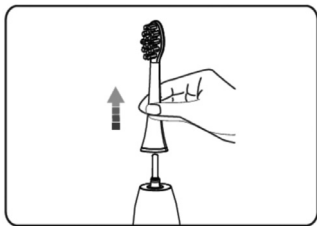- Remove and rinse the toothbrush head thoroughly.
- To clean the shaft and body of the appliance, run a little water over them. To do this, hold the appliance body so that the shaft is pointing down.
- Do not wash with hot water above 40°C. Do not use caustic or abrasive products.
- Wipe the appliance body with a towel or cloth.
- Replace the toothbrush head if the bristles start to bend outwards.

## **STORAGE**

Before storage clean the appliance according to the «Cleaning and care» section. Avoid exposure to extreme temperatures, direct sunlight and excessive humidity. Do not store in a dusty environment. Keep the appliance clean.

## **TROUBLESHOOTING**

<b>Problem</b>	<b>Possible reason</b>	<b>Possible solution</b>
The appliance does not turn on	Low battery level	Charge the appliance
The toothbrush head cannot be installed on the appliance body	Foreign objects in the toothbrush head opening	Remove foreign objects
	The shaft is dirty	Clean the shaft from dirt
	The toothbrush head is deformed	Replace the toothbrush head
The appliance operates for a short period of time	Insufficient battery charge	Charge the appliance fully
	Battery life has expired	Contact an authorized service center
Low vibration power	The toothbrush head is installed incorrectly or is dirty	Clean the inside of the toothbrush head and the shaft on the appliance body from dirt, and reinstall the toothbrush head
	The toothbrush head is pressed too hard against the teeth	When brushing lightly press the toothbrush head against the surface of the teeth
	Low battery level	Charge the appliance

**If the above actions did not help to eliminate the problem do not try to repair the appliance by yourself, contact an authorized service center for help.**

## **APPLIANCE PRODUCTION DATE**

The production date is indicated on the package (in the serial number) after «S/N», in the format: MMYYYYXXXXXXXXXX, where MM is the month of production, YYYY is the year of production. Do not throw away the packaging before the warranty period has expired.

## **CORRECT DISPOSAL OF THE APPLIANCE**

This marking indicates that this appliance should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of resources. Dispose of the appliance through an authorized disposal company or an appropriate local authority.

A recyclable lithium-ion battery powers this appliance. When the appliance reaches the end of its useful life, remove the battery from the appliance before disposing of it. Dispose of the battery in accordance with local disposal requirements.

Dispose of the packaging materials environmentally. The packaging is made of recyclable materials which should be disposed of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning: 1–7: plastics / 20–22: paper and fiberboard / 80–98: composite materials.

If in doubt, contact waste disposal company. Local community or municipal government can provide information on how to dispose of the appliance properly.

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## Warranty conditions:

---

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 3 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid.

Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by Authorized Service Center or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card. This warranty does not cover installation and uninstallation of the appliance. All costs associated with the installation and uninstallation of the appliance (including warranty repairs cases) are borne by the customer.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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# WARRANTY CARD

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this product.

**Keep this card during the warranty period. When purchasing a product ask for the warranty card to be completely filled out.**

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. \*
2. The service life of household appliances is 3 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
  - misuse and non-consumer use;
  - mechanical damage;
  - damage caused by the ingress of foreign objects, substances, liquids, insects;
  - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
  - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
  - when operating equipment in the power supply network with a missing single ground loop;
  - in case of violation of seals installed on the product;
  - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories as well as blades, filters, shelves, drawers, stands and other containers for storing products, etc.
8. The warranty period is 12 months from the date of sale.

\* Tear-off maintenance tickets are provided by an authorized service center.

**The completeness of the product is checked. I have read the terms of the warranty service, no complaints.**

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